

OLYMPUS

Your Vision, Our Future

ETS PLUS

Secure Endoscope Transportation with Hygienic Focus



IMPROVE THE EFFICIENCY OF YOUR DAILY WORK

The endoscope transport system “ETS Plus” is an advanced and dedicated product to simultaneously transport up to six flexible medical endoscopes between the examination room, the reprocessing area and the storage area.

Endoscope Transport – User Friendly, Simple and Safe

The ETS Plus is specifically designed to provide easy and swift transport of all flexible medical endoscopes.

1 Ergonomic Height

Loading and unloading at a user friendly working height is supported by the top shelf of the ETS Plus, allowing safe packing and unpacking of the endoscopes and avoiding the necessity for additional worktops.

2 Additional Space

The ETS Plus comes with six trays that are placed on tilted slide rails, tightly securing the valuable content during transport.

3 Tray Protection

The stainless steel frame perfectly encloses and protects all the trays.

4 Improved Handling

Comfortable transport of the endoscopes is assured by the ETS Plus handle and the base plate cut-out of the trolley.


5 Easy Movement

Four freely rotatable wheels, some fitted with brakes, enable full and easy control when moving through different departments.



Hygienic Safety

The separately packed foil kit consists of two robust foils: an inner liner to avoid contamination of the tray and a cover to protect against dust and recontamination. The cover is printed on both sides, indicating the status of the endoscope in the tray by a red/green color code and a symbol. The transparency of the cover allows easy visual identification of the endoscope placed in the tray.

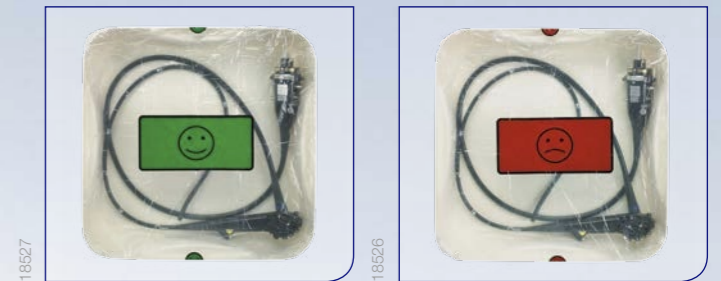

Quick Foil Placement

Placing an endoscope with inner liner and cover from a foil kit is an easy three-step routine. Both foils are tightly fixed to the tray. The single-use application further reduces the risk of contamination during transport.



Easy and Effective Cleaning

The plain surfaces of the ETS Plus trays offer compatibility for all types of flexible medical endoscopes while being easy to clean and disinfect. The trays have no overhanging edges, and all four grips face inwards, reducing places for bacteria to grow.



Specifications and Technical Data

Specifications	Number of instrument trays	6
	Number of grips per tray	4
	Wheels	4 (360°C freely rotatable; 2 with brake)
	Application	In-house transport of medical flexible endoscopes
Dimensions	ETS Plus – Transfer trolley	W x D x H: 600 x 470 x 1050 mm Weight: 30 kg
	ETS Plus – Instrument tray	W x D x H: 520 x 520 x 100 mm Weight: < 1 kg
	ETS Plus - Wall dispenser	W x D x H: 150 x 130 x 420 mm Weight: < 1 kg
	Working height	1050 mm
Disinfection	Transfer trolley	Wipe disinfection (use common hospital surface disinfectant agents)
	Instrument tray	Automated thermal disinfection (min. A ₀ 600) Thermostability: ≤ 93°C
Foil Kit	Contents	1x Liner foil 1x Cover foil
	Single use	Yes
	Sterile	No
Standards	Medical device classification	No

Specifications, design, and accessories are subject to change without any notice or obligation on the part of the manufacturer.